

# Ficolin-1 ELISA Kit (Human) (OKBB00565)

## **Instruction for Use**

For the quantitative measurement of Human Ficolin-1 in cell culture supernatants, serum and plasma (heparin).

This product is intended for research use only.



## Table of Contents

| 1.  | Background                      | . 2 |
|-----|---------------------------------|-----|
| 2.  | Assay Summary                   | 2   |
| 3.  | Precautions                     | . 3 |
| 4.  | Storage and Stability           | . 3 |
| 5.  | Kit Components                  | . 3 |
| 6.  | Required Materials Not Supplied | . 3 |
| 7.  | Technical Application Tips      | . 4 |
| 8.  | Reagent Preparation             | . 4 |
| 9.  | Sample Preparation              | . 6 |
| 10. | Assay Procedure                 | . 7 |
| 11. | Calculation of Results          | . 8 |
| 12. | Typical Expected Data           | . 8 |
| 13. | Technical Resources             | 10  |



## 1. Background


#### **Principle**

Aviva Systems Biology Ficolin-1 ELISA Kit (Human) (OKBB00565) is based on standard sandwich enzymelinked immuno-sorbent assay technology. A mouse monoclonal antibody specific for Ficolin-1 has been precoated onto 96-wellplate (12 x 8 Well Strips). Standards (NSO; A30-A326) and test samples are added to the wells and incubated. After washing, a biotinylated polyclonal goat detector antibody specific for Ficolin-1 is added, incubated and followed by washing. Avidin-Biotin-Peroxidase Complex is then added, incubated and unbound conjugate is washed away. An enzymatic reaction is visualized through the addition of TMB substrate which is catalyzed by HRP to produce a blue color product that changes yellow after adding acidic stop solution. The density of yellow coloration read by absorbance at 450 nm and is quantitatively proportional to the amount of sample Human Ficolin-1 captured in well.

#### **Background**

Ficolin-1, and also commonly termed M-ficolin is a protein that in humans is encoded by the FCN1 gene. It is mapped to 9q34.3. The ficolin family of proteins are characterized by the presence of a leader peptide, a short N-terminal segment, followed by a collagen-like region, and a C-terminal fibrinogen-like domain. Ficolins have a crucial role in defense against pneumococcal infection through the lectin complement pathway. Ficolin 1 is predominantly expressed in the peripheral blood leukocytes, and has been postulated to function as a plasma protein with elastin-binding activity.

## 2. Assay Summary

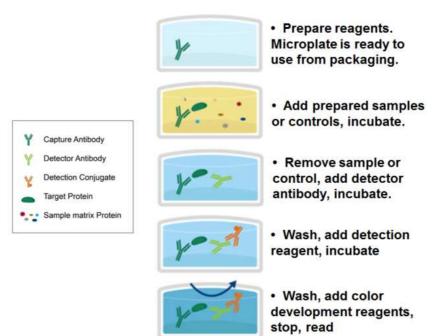



## 3. Precautions

- Read instructions fully prior to beginning use of the assay kit.
- Any deviations or modifications from the described method or use of other reagents could result in a reduction of performance.
- Reduce exposure to potentially harmful substances by wearing personal protective lab equipment including lab coats, gloves and glasses.
- For information on hazardous substances included in the kit please refer to the Material Safety Data Sheet (MSDS).

## 4. Storage and Stability

• Upon receipt store kit at 4°C for 6 months or -20°C for 12 months. Avoid multiple freeze/thaw cycles.

## 5. Kit Components

• The following reagents are the provided contents of the kit.

| Description   | Quantity              |  |  |
|---|-----------------------|--|--|
| 96-Well plate Pre-coated with Anti-Human Ficolin-1 Antibody | 1 (12 x 8 Well Strip) |  |  |
| Lyophilized Recombinant Human Ficolin-1 Standard            | 50 ng/tube × 2        |  |  |
| 100X Biotinylated Anti-Human Ficolin-1 Antibody             | 130 µL                |  |  |
| 100X Avidin-Biotin-Peroxidase Complex (ABC)                 | 130 µL                |  |  |
| Sample Diluent Buffer                                       | 30 mL                 |  |  |
| Antibody Diluent Buffer                                     | 12 mL                 |  |  |
| ABC Diluent Buffer  | 12 mL                 |  |  |
| TMB Color Developing Agent                                  | 10 mL                 |  |  |
| TMB Stop Solution   | 10 mL                 |  |  |
| 10X Wash Buffer   | 30 mL                 |  |  |

## 6. Required Materials Not Supplied

- Microplate reader capable of reading absorbance at 450 nm.
- Automated plate washer (optional).
- Pipettes capable of precisely dispensing 0.5 µL through 1 mL volumes of aqueous solutions.
- Pipettes or volumetric glassware capable of precisely measuring 1 mL through 100 mL of aqueous solutions.
- New, clean tubes and/or micro-centrifuge tubes for the preparation of standards or samples.
- Absorbent paper or paper toweling.
- Distilled or deionized ultrapure water.



## 7. Technical Application Tips

- Do not mix or substitute components from other kits.
- To ensure the validity of experimental operation, it is recommended that pilot experiments using standards and a small selection of sample dilutions to ensure optimal dilution range for quantitation.
- Samples exhibiting OD measurements higher than the highest standard should be diluted further in the appropriate sample dilution buffers.

• Inspect all reagents prior to use. Components should contain no particulates or cloudiness and should be colorless.

- Prior to using the kit, briefly spin component tubes to collect all reagent at the bottom.
- Replicate wells are recommended for standards and samples.
- Cover microplate while incubating to prevent evaporation.
- Do not allow the microplate wells dry at any point during the assay procedure.
- Do not reuse tips or tube to prevent cross contamination.
- Avoid causing bubbles or foaming when pipetting, mixing or reconstituting.
- Completely remove of all liquids when washing to prevent cross contamination.
- Prepare reagents immediately prior to use and do not store, with the exception of the top standard.
- Equilibrate all materials to ambient room temperature prior to use (standards exception).

• For optimal results in inter- intra-assay consistency, equilibrate all materials to 37°C prior to performing assay (standards exception) and perform all incubations at 37°C.

 $\bullet$  Pipetting less than 1  $\mu L$  is not recommended for optimal assay accuracy.

• Once the procedure has been started, all steps should be completed without interruption. Ensure that all reagents, materials and devices are ready at the appropriate time.

• Incubation times will affect results. All wells should be handled in the same sequential order and time intervals for optimal results.

• Samples containing precipitates or fibrin strands or which are hemolytic of lipemic might cause inaccurate results due to interfering factors.

## 8. Reagent Preparation

• Equilibrate all materials to room temperature prior to use and use prepare immediately prior to use.

#### 8.1 1X Biotinylated Anti-Human Ficolin-1 Antibody

- 8.1.1 Prepare the **1X Biotinylated Anti-Human Ficolin-1 Antibody** immediately prior to use by diluting the **100X Biotinylated Anti-Human Ficolin-1 Antibody** 1:100 with Antibody Diluent Buffer.
- 8.1.2 For each well to be used in the experiment prepare 120 μL by adding 1.2 μL of 100X
  Biotinylated Anti-Human Ficolin-1 Antibody to 118.8 μL Antibody Diluent Buffer.
- 8.1.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure.

#### 8.2 1X Avidin-Biotin-Peroxidase Complex (ABC)


- 8.2.1 Prepare the **1X Avidin-Biotin-Peroxidase Complex (ABC)** immediately prior to use by diluting the **100X Avidin-Biotin-Peroxidase Complex (ABC)** 1:100 with **ABC Dilution Buffer**.
- 8.2.2 For each well to be used in the experiment prepare 120 μL, by adding 1.2 μL of **100X Avidin-Biotin-Peroxidase Complex (ABC)** to 118.8 μL **ABC Dilution Buffer**.
- 8.2.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure.



#### 8.3 Ficolin-1 Assay Standards

- 8.3.1 Prepare the Ficolin-1 standards no greater than 2 hours prior to performing experiment. Standards should be held on ice until use in the experiment.
- 8.3.2 Reconstitute one of the provided 10 ng Lyophilized Recombinant Human Ficolin-1 Standard. Use one for each experiment. Prepare a stock 10,000 pg/mL Human Ficolin-1 Standard by reconstituting one tube of Lyophilized Recombinant Human Ficolin-1 Standard as follows:
  8.3.2.1 Gently spin or tap the vial to collect all material at the bottom.
  - 0.3.2.1 Gentily spin of tap the vial to collect all material at the r
  - 8.3.2.2 Add 1 mL of Sample Diluent Buffer to the vial.
  - 8.3.2.3 Seal then mix gently and thoroughly.
  - 8.3.2.4 Leave the vial to sit at ambient temperature for 10 minutes.
- 8.3.3 Prepare a set of seven serially diluted standards as follows:
  - 8.3.3.1 Label tubes with numbers 2 8.
    - 8.3.3.2 Add 300  $\mu$ L of **Sample Diluent Buffer** to Tube #'s 2 8.
    - 8.3.3.3 Use the reconstituted 50 ng/mL Human Ficolin-1 as Standard #1.
    - 8.3.3.4 Prepare Standard #2 by adding 300 μL of 50 ng/mL Standard #1 from Tube #1 to Tube #2. Mix gently and thoroughly.
    - 8.3.3.5 Prepare **Standard #3** by adding 300 μL of **Standard #2** from Tube #2 to Tube #3. Mix gently and thoroughly.
    - 8.3.3.6 Prepare further serial dilutions through Tube #7. Reference the table below as a guide for serial dilution scheme.
    - 8.3.3.7 Tube #8 is a blank standard (only **Sample Diluent Buffer**), which should be included with every experiment.

| Standard<br>Number<br>(Tube) | Sample To Dilute                        | Volume<br>Standard<br>(µL) | Volume Sample<br>Diluent Buffer<br>(µL) | Total Volume<br>(μL) | Final Concentration |
|------------------------------|---|----------------------------|---|----------------------|---------------------|
| 1                            | 50 ng/mL of Human<br>Ficolin-1 standard | NA                         | NA                                      | 1,000                | 50 ng/mL            |
| 2                            | Tube #1                                 | 300                        | 300                                     | 600                  | 25 ng/mL            |
| 3                            | Tube #2                                 | 300                        | 300                                     | 600                  | 12.5 ng/mL          |
| 4                            | Tube #3                                 | 300                        | 300                                     | 600                  | 6.25 ng/mL          |
| 5                            | Tube #4                                 | 300                        | 300                                     | 600                  | 3.12 ng/mL          |
| 6                            | Tube #5                                 | 300                        | 300                                     | 600                  | 1.56 ng/mL          |
| 7                            | Tube #6                                 | 300                        | 300                                     | 600                  | 0.78 ng/mL          |
| 8                            | NA                                      | 0                          | 300                                     | 300                  | 0.0 (Blank)         |

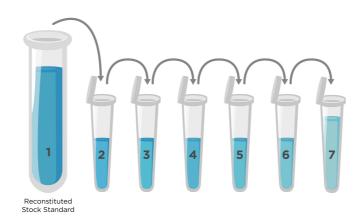



#### 8.4 1X Wash Buffer

- 8.4.1 Add 270 mL of ultra-pure water to a clean > 500 mL bottle or other vessel.
- 8.4.2 Add the entire 30 mL contents of the **10X Wash Buffer** bottle to the water.
- 8.4.3 Seal and mix gently by inversion. Avoid foaming or bubbles.
- 8.4.4 Store the **1X Wash Buffer** at room temperature until ready to use in the procedure. Store the prepared **1X Wash Buffer** at 4°C for no longer than 1 week. Do not freeze.

#### 8.5 Microplate Preparation

- Micro-plates are provided ready to use and do not require rinsing or blocking.
- Unused well strips should be returned to the original packaging, sealed and stored at 4°C.
- Equilibrate microplates to ambient temperatures prior to opening to reduce potential condensation.

## 9. Sample Preparation

#### 9.1 Sample Preparation and Storage

- Store samples to be assayed at 2-8°C for 24 hours prior being assayed.
- For long term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.
- Clear samples by centrifugation as follows:

• Cell culture supernatant - Remove particulates by centrifugation, assay immediately or aliquot and store samples at -20°C.

• Serum - Allow the serum to clot in a serum separator tube (about 4 hours) at room temperature. Centrifuge at approximately 1,000 x g for 15 min. Analyze the serum immediately or aliquot and store samples at -20°C.

• **Plasma** - Collect plasma using heparin as an anticoagulant. Centrifuge for 15 min at 1,500 x g within 30 min of collection. Assay immediately or aliquot and store samples at -20°C.

#### 9.2 Sample Dilution

Target protein concentration must be estimated and appropriate sample dilution selected such that the final target protein concentration falls near the middle of the assay linear dynamic range.

- $\bullet$  Prepare 150  $\mu L$  sample for each replicate to be assayed.
- Dilute samples with Sample Diluent Buffer.
- Mix diluted samples gently and thoroughly.
- Pipetting less than 2 µL is not recommended for optimal assay accuracy.

• Refer to the following table for recommended sample dilution guidelines based on the dynamic range of this kit:

| Estimated Sample T     | arget Concentration | Dilution Level     | Sample Volume<br>For Two<br>Replicates | Sample Diluent Buffer<br>For Two Replicates |  |
|------------------------|---------------------|--------------------|--|---|--|
| High Concentration     | 500-5,000 ng/mL     | 1:100              | 3 µL                                   | 297 µL                                      |  |
| Medium Concentration   | 50-500 ng/mL        | 1:10               | 30 µL                                  | 270 µL                                      |  |
| Low Concentration      | 0.78-50 ng/mL       | 1:2                | 150 µL                                 | 150 μL                                      |  |
| Very Low Concentration | ≤0.78 ng/mL         | 1:2 or No Dilution | -                                      | -   |  |



## **10. Assay Procedure**

- Equilibrate all reagents and materials to ambient room temperature prior to use in the procedure.
- To control for small potential variations in micro well-plate and day to day ambient temperature fluctuations, equilibrate all reagents prior to use and perform all incubation steps at to 37°C for optimal consistency and reproducibility.
- **10.1** Add 100 µL of serially titrated standards, diluted samples or blank into wells of the pre-coated well plate. At least two replicates of each standard, sample or blank is recommended.
- **10.2** Cover the plate with the well plate lid and incubate for 90 minutes.
- **10.3** Remove the plate lid and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
- **10.4** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- **10.5** Add 100 µL of prepared **1X Biotinylated Anti-Human Ficolin-1 Antibody** to each well.
- **10.6** Cover with the well-plate lid and Incubate for 60 minutes.
- **10.7** Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
- **10.8** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 10.9 Wash plate 3 times with 1X Wash Buffer as follows:
  - 10.9.1 Add 300 µL of 1X Wash Buffer to each assay well.
  - 10.9.2 Incubate for 1 minute.
  - 10.9.3 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
  - 10.9.4 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
  - 10.9.5 Repeat steps 10.9.1 through 10.9.4 two more times
- **10.10** Add 100 μL of prepared **1X Avidin-Biotin-Peroxidase Complex (ABC)** into each well and incubate for 30 minutes.
- **10.11** Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
- **10.12** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 10.13 Wash plate 5 times with 1X Wash Buffer as follows:
  - 10.13.1 Add 300  $\mu L$  of 1X Wash Buffer to each assay well.
  - 10.13.2 Incubate for 1 minute.
  - 10.13.3 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
  - 10.13.4 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
  - 10.13.5 Repeat steps 10.13.1 through 10.13.4 four more times.
- **10.14** Add 90 μL of **TMB Color Developing Agent** to each well and incubate in the dark for 25 30 minutes. (NOTE: optimal incubation time must be determined by the user. Optimal development can be visualized by blue shading in the top four standard wells, while the remaining standards are still clear.)
- 10.15 Add 100  $\mu L$  of TMP Stop Solution to each well. Well color should change to yellow immediately.
- **10.16** Read the O.D. absorbance at 450 nm with a standard microplate reader within 30 minutes of stopping the reaction in step 10.15.



## 11. Calculation of Results

For analysis of the assay results, calculate the Relative OD<sub>450</sub> for each test or standard well as follows:

(Relative  $OD_{450}$ ) = (Well  $OD_{450}$ )–(Mean Blank Well  $OD_{450}$ )

The standard curve is generated by plotting the mean replicate **Relative OD**<sub>450</sub> of each standard serial dilution point vs. the respective standard concentration. The Human Ficolin-1 concentration contained in the samples can be interpolated by using linear regression of each mean sample **Relative OD**<sub>450</sub> against the standard curve. This is best achieved using curve fitting software.


**Note:** if the samples measured were diluted, multiply the derived mean sample concentration by the dilution factor for a final sample concentration.

## 12. Typical Expected Data

12.1 <u>Typical absorbance values</u>. Expected absorbance for standards when TMB incubation is performed for 20 minutes at 37°C and measured at OD<sub>450</sub>.

| Standard<br>Number       | 8     | 7     | 6     | 5     | 4    | 3     | 2     | 1     |
|--------------------------|-------|-------|-------|-------|------|-------|-------|-------|
| Concentration<br>(pg/mL) | 0     | 0.78  | 1.56  | 3     | 6    | 13    | 25    | 50    |
| OD <sub>450</sub>        | 0.066 | 0.187 | 0.323 | 0.557 | 0.93 | 1.423 | 1.967 | 2.447 |

**12.2** <u>Typical standard curve</u>. This standard curve is for demonstration purposes only. An assay specific standard curve should be performed with each assay.

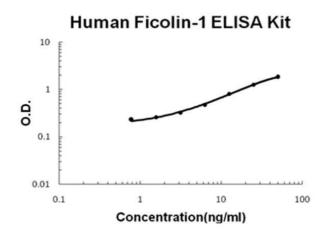



## 12.3 General Specifications

| General Specificiations |  |  |  |  |  |  |
|-------------------------|--|--|--|--|--|--|
| Range                   | 0.78 ng/mL-50 ng/mL  |  |  |  |  |  |
| Sensitivity             | < 20 pg/mL (Derived by linear regression of $OD_{450}$ of the Mean Blank + 2xSD) |  |  |  |  |  |
| Specificity             | Natural and recombinant Human Ficolin-1<br>UniProt ID: O00602                    |  |  |  |  |  |
| Cross-Reactivity        | No detectable cross-reactivity with other relevant proteins                      |  |  |  |  |  |
| Recovery                | -  |  |  |  |  |  |
| Linearity               | -  |  |  |  |  |  |

## 12.4 <u>Reproducibility</u>

|  | Intra-Assay |       |       | Inter-Assay |      |      |
|--|-------------|-------|-------|-------------|------|------|
| Sample ID                              | 1           | 2     | 3     | 1           | 2    | 3    |
| n =                                    | 16          | 16    | 16    | 24          | 24   | 24   |
| Mean Measured<br>Concentration (pg/mL) | 4.83        | 13.2  | 25.1  | 5.24        | 12.9 | 25.9 |
| Standard Deviation (pg/mL)             | 0.188       | 0.581 | 1.431 | 0.393       | 0.89 | 1.89 |
| Consistency (%CV)                      | 3.9         | 4.4   | 5.7   | 7.5         | 6.9  | 7.3  |



## 13. Technical Resources

#### 13.1 References

- 13.1.1 Endo, Y., Takahashi, M., Iwaki, D., Ishida, Y., Nakazawa, N., Kodama, T., Matsuzaka, T., Kanno, K., Liu, Y., Tsuchiya, K., Kawamura, I., Ikawa, M., Waguri, S., Wada, I., Matsushita, M., Schwaeble, W. J., Fujita, T. Mice deficient in ficolin, a lectin component pathway recognition molecule, are susceptible to Streptococcus pneumoniae infection.
- 13.1.2 "Entrez Gene: FCN1 ficolin (collagen/fibrinogen domain containing) 1" J. Immun. 189: 5860-5866, 2012.
- 13.1.3 Lu J, Tay PN, Kon OL, Reid KB (Mar 1996). "Human ficolin: cDNA cloning, demonstration of peripheral blood leucocytes as the major site of synthesis and assignment of the gene to chromosome 9".Biochem J 313 (2): 473–8.
- **13.2** Technical Support

USA

Aviva Systems Biology, Corp. 5754 Pacific Center Blvd, Suite 201 San Diego, CA 92121

Phone: 858-552-6979 Toll Free: 888-880-0001 Fax: 858-552-6975

Technical support: techsupport@avivasysbio.com

<u>China</u>

Beijing AVIVA Systems Biology 6th Floor, B Building, Kaichi Tower #A-2 Jinfu Road. Daxing Industrial Development Zone Beijing, 102600, CHINA

Phone: (86)10-60214720 Fax: (86)10-60214722 E-mail: support@avivasysbio.com.cn

中国地址:北京大兴工业开发区金辅路甲2号凯驰大厦B座6层(102600)

电话: 010-60214720/21 传真: 010-60214722

产品售前咨询及销售: sales@avivasysbio.com.cn

售后及技术支持: support@avivasysbio.com.cn