



**Ciprofloxacin ELISA Kit
(OKAO00115)
Lot# KH0689**

Instructions for Use

For the quantitative detection of Ciprofloxacin (CPFX) concentration in tissues, honey, milk, milk powder and egg.

Variation between lots can occur. Refer to the manual provided with the kit.

This product is intended for research use only.

Table of Contents

| | |
|--|----|
| 1. Background | 2 |
| 2. Storage and Stability..... | 3 |
| 3. Kit Components | 3 |
| 4. Required Materials Not Supplied | 3 |
| 5. Precautions | 4 |
| 6. Technical Application Tips | 4 |
| 7. Reagent Preparation..... | 5 |
| 8. Sample Preparation | 6 |
| 9. Assay Procedure | 7 |
| 10. Calculation of Results..... | 8 |
| 11. Typical Expected Data | 8 |
| 12. Technical Resources | 10 |

1. Background

Principle

Aviva Systems Biology Ciprofloxacin ELISA Kit (OKAO00115) is based on standard competitive inhibition enzyme-linked immuno-sorbent assay technology. Ciprofloxacin has been pre-coated onto a 96-wellplate (12 x 8 Well Strips). Standards or test samples are added to the wells along with an anti-Ciprofloxacin antibody and an HRP Conjugated antibody, then incubated. Wells are washed and detection substrates are added and incubated. An enzymatic reaction is produced through the addition of TMB which is catalyzed by HRP generating a blue color product that changes yellow after adding acidic stop solution. The density of yellow coloration read by absorbance at 450 nm which is quantitatively proportional to the amount of detector antibody bound in the well and inversely proportional to the amount of Ciprofloxacin in the sample.

General Specifications

| General Specifications | |
|------------------------|--|
| Range | 0.1 – 8.1 ppb |
| LOD | < 0.1 ppb (Derived by linear regression of OD ₄₅₀ of the Mean Blank + 2xSD) |
| Specificity | Ciprofloxacin <u>PubChem#</u> : 2764 <u>Chemical Formula</u> : C ₁₇ H ₁₈ FN ₃ O ₃ <u>Alias</u> : CPFX; Ciprofloxacin; Ciprobay; 85721-33-1; Ciproxan; Ciprofloxacin; Ciprinol |

2. Storage and Stability

- Upon receipt store kit at 4°C. Do not use past expiration date.

3. Kit Components

- The following reagents are the provided contents of the kit.

| Description | Quantity | Storage Conditions |
|---------------------------------|-------------------------------|---|
| Ciprofloxacin Microplate | 96 Wells (12 x 8 Well strips) | Store at 4°C Do not use past expiration date |
| Ciprofloxacin Standards | 6 x 1 mL | |
| Ciprofloxacin Detector Antibody | 1 x 6 mL | |
| HRP Conjugate | 1 x 6 mL | |
| Detection Reagent A | 1 x 6 mL | |
| Detection Reagent B | 1 x 6 mL | |
| 5X Sample Prep Buffer | 1 x 50 mL | |
| 20X Wash Buffer | 1 x 20 mL | |
| Stop Solution | 1 x 6 mL | |

4. Required Materials Not Supplied

- Microplate reader capable of reading absorbance at 450 nm.
- Automated plate washer (optional).
- Pipettes capable of precisely dispensing 0.5 µL through 1 mL volumes of aqueous solutions.
- Pipettes or volumetric glassware capable of precisely measuring 1 mL through 100 mL of aqueous solutions.
- New, clean tubes and/or micro-centrifuge tubes for the preparation of standards or samples.
- Absorbent paper or paper toweling.
- Distilled or deionized ultrapure water.
- Heater water bath
- Centrifuge
- Reagents: dichloromethane, n-Hexane, HCl (approximately 36.5%), Acetonitrile (CH₃CN) 99.9%.

5. Precautions

- Read instructions fully prior to beginning use of the assay kit.
- Any deviations or modifications from the described method or use of other reagents could result in a reduction of performance.
- Reduce exposure to potentially harmful substances by wearing personal protective lab equipment including lab coats, gloves and glasses.
- For information on hazardous substances included in the kit please refer to the Material Safety Data Sheet (MSDS).
- Kit cannot be used beyond the expiration date on the label.

6. Technical Application Tips

- Do not mix or substitute components from other kits.
- To ensure the validity of experimental operation, it is recommended that pilot experiments using standards and a small selection of sample dilutions to ensure optimal dilution range for quantitation.
- Samples exhibiting OD measurements higher than the highest standard should be diluted further in the appropriate sample dilution buffers.
- Prior to using the kit, briefly spin component tubes to collect all reagents at the bottom.
- Replicate wells are recommended for standards and samples.
- Cover microplate while incubating to prevent evaporation.
- Do not allow the microplate wells dry at any point during the assay procedure.
- Do not reuse tips or tube to prevent cross contamination.
- Avoid causing bubbles or foaming when pipetting, mixing or reconstituting.
- Completely remove of all liquids when washing to prevent cross contamination.
- Prepare reagents immediately prior to use and do not store, with the exception of the top standard.
- Equilibrate all materials to ambient room temperature prior to use (standards exception).
- For optimal results for inter- and intra-assay consistency, equilibrate all materials to room temperature prior to performing assay (standards exception) and perform all incubations at 37°C.
- Pipetting less than 1 μ L is not recommended for optimal assay accuracy.
- Once the procedure has been started, all steps should be completed without interruption. Ensure that all reagents, materials and devices are ready at the appropriate time.
- Incubation times will affect results. All wells should be handled in the same sequential order and time intervals for optimal results.
- Samples containing precipitates, fibrin strands or bilirubin, or are hemolytic or lipemic might cause inaccurate results due to interfering factors.
- TMB Substrate is easily contaminated and should be colorless or light blue until added to plate. Handle carefully and protect from light.
- To minimize influences on the assay performance, operation procedures and lab conditions, especially room temperature, air humidity, incubator temperature should be strictly controlled. It is also strongly suggested that the whole assay is performed by the same operator from the beginning to the end.

7. Reagent Preparation

- Equilibrate all materials to room temperature prior to use and use prepare immediately prior to use.
- The following reagents are provided at ready to use concentrations and require no preparation:
 - Ciprofloxacin Standards
 - Ciprofloxacin Detector Antibody
 - HRP Conjugate
 - Detection Reagent A
 - Detection Reagent B

7.1 Standards

Standards are provided at ready to use concentrations. The concentrations are as follows:

| Tube | Standard 5 | Standard 4 | Standard 3 | Standard 2 | Standard 1 | Standard 0 |
|------|---------------|---------------|---------------|---------------|---------------|---------------|
| ppb | 8.1 | 2.7 | 0.9 | 0.3 | 0.1 | 0 |

7.2 1X Wash Buffer

- 7.2.1 Add the 20 mL of **20X Wash Buffer** to 380 mL of ultrapure water.
- 7.2.2 Mix gently and thoroughly. Store at 4°C for one week.

7.3 1X Sample Prep Buffer

- 7.3.1 Add the 50 mL of **5X Sample Prep Buffer** to 200 mL of ultrapure water.
- 7.3.2 Mix gently and thoroughly. Store at 4°C for one week.

7.4 0.1 M HCl

- 7.4.1 Prepare 100 mL of fresh **0.1 M HCl** immediately prior to use.
- 7.4.2 Add 860 µL HCl (36%) to 100 mL ultrapure water.
- 7.4.3 Mix gently and thoroughly.
- 7.4.4 Do not store for future use.

7.5 Acetonitrile-HCL

- 7.5.1 Prepare an **Acetonitrile-HCL** solution by combining Acetonitrile and 0.1M HCL in a ratio of 4:1 (Acetonitrile : 0.1M HCL)

7.6 N-Hexane-Dichloromethane

- 7.6.1 Prepare an **N-Hexane-Dichloromethane** solution by combining N-Hexane and Dichloromethane in a ratio of 1:4 (N-Hexane : Dichloromethane).

7.7 Microplate Preparation

- Micro-plates are provided ready to use and do not require rinsing or blocking.
- Unused well strips should be returned to the original packaging, sealed and stored at 4°C.
- Equilibrate microplates to ambient temperatures prior to opening to reduce potential condensation.

8. Sample Preparation

8.1 Sample Preparation and Storage

- Mix diluted samples gently and thoroughly.
- Pipetting less than 2 µL is not recommended for optimal assay accuracy.
- Untreated samples should be stored at -20°C.
- Prepared samples can be stable for 1 week at 4°C.

8.1.1 **Tissue (chicken, liver, shrimp, fish, etc.) (Dilution 1:2)**

- 8.1.1.1 Homogenize the sample at 10,000 rpm for 1 minute.
- 8.1.1.2 Weigh 2.0 ± 0.05 g of the homogenized sample into a 50 mL centrifugal tube, add 8 mL **Acetonitrile-HCL**.
- 8.1.1.3 Shake for vigorously for 5 minutes.
- 8.1.1.4 Centrifuge at above 4,000 rpm at room temperature (20-25 °C) for 10 minutes.
- 8.1.1.5 Transfer 2 mL of supernatant into a new centrifugal tube, dry completely with a gentle stream of nitrogen or rotary evaporation at 56 °C.
- 8.1.1.6 Add 1 mL **N-hexane** to dry residue, shake for 2 minutes.
- 8.1.1.7 Add 1 mL of the diluted **1X Sample Prep Buffer**, shake strongly for 30 seconds.
- 8.1.1.8 Centrifuge at 4,000 rpm at room temperature for 5 minutes, remove the upper layer.
- 8.1.1.9 Assay the lower layer sample neat (do not dilute).

8.1.2 **Honey (Dilution 1:2)**

- 8.1.2.1 Weight 1.0 ± 0.05 g honey into 50 mL centrifugal tube, then add 6 mL **Acetonitrile-HCL**. Shake properly for 5 minutes to dissolve completely.
- 8.1.2.2 Add 3 mL of **1X Sample Prep Buffer**, then add 11 mL **Dichloromethane**.
- 8.1.2.3 Shake for 5 minutes then centrifuge greater than 4,000 rpm at room temperature (20-25°C) for 5 minutes.
- 8.1.2.4 Remove the upper layer.
- 8.1.2.5 Transfer 8 mL of the lower layer to a fresh 50 mL centrifuge tube.
- 8.1.2.6 Dry completely with nitrogen by rotary evaporation at 56 °C.
- 8.1.2.7 Add 1 mL of **1X Sample Prep Buffer** and shake for 2 minutes.
- 8.1.2.8 Add 1 mL of **N-hexane** and mix for 30 seconds.
- 8.1.2.9 Centrifuge greater than 3,000 rpm at room temperature (20-25°C) for 5 minutes.
- 8.1.2.10 Discard the upper layer. Assay 50 µL of the lower layer samples.

8.1.3 **Egg (Diluted 1:30)**

- 8.1.3.1 Homogenize the sample.
- 8.1.3.2 Weigh 1.0 ± 0.02 g of the homogenized sample into 10 mL centrifugal tube.
- 8.1.3.3 Add 5 mL of deionized water, shake for 1 minute to dissolve completely.
- 8.1.3.4 Transfer 100 µL sample solution into a 2 mL centrifuge tube.
- 8.1.3.5 Add 400 µL of **1X Sample Prep Buffer**, shake for 1 minute.
- 8.1.3.6 Take 50 µL for further analysis.

8.1.4 **Milk (Diluted 1:20)**

- 8.1.4.1 Take 5 mL sample into 50 mL centrifugal tube.
- 8.1.4.2 Centrifuge at above 4,000 rpm at 15°C for 10 minutes. Remove fat (upper layer).
- 8.1.4.3 Transfer 25 µL milk removed fat into 2 mL centrifuge tube.

8.1.4.4 Add 475 μ L of the **1X Sample Prep Buffer**, shake for 1 minute to dissolve completely.

8.1.4.5 Take 50 μ L for further analysis.

8.1.5 Milk Powder (Diluted 1:50)

8.1.5.1 Weigh 0.5 ± 0.02 g milk powder into 10 mL centrifuge tube.

8.1.5.2 Add 5 mL of deionized water, shake for 1 minute to dissolve completely.

8.1.5.3 Transfer 100 μ L of sample solution into 2 mL centrifuge tube.

8.1.5.4 Add 400 μ L of **1X Sample Prep Buffer**, shake for 1 minute.

8.1.5.5 Take 50 μ L for further analysis.

9. Assay Procedure

- Equilibrate all reagents and materials to ambient room temperature prior to use in the procedure.
- Temperature regulation for all procedures and incubations to 25°C is recommended for optimal reproducibility.

- 9.1 Determine the required number of wells and return any remaining unused wells and desiccant to the pouch.
- 9.2 Add 50 μ L of standards or samples into wells of the **Ciprofloxacin Microplate**. At least two replicates of each standard, sample or blank is recommended.
- 9.3 Immediately add 50 μ L of the **Ciprofloxacin Detector Antibody** and 50 μ L of the **HRP Conjugate** to each well.
- 9.4 Cover the plate with the well plate sealer, gently mix and incubate for 45 minutes.
- 9.5 Remove the plate sealer and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
- 9.6 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 9.7 Wash plate 5 times with **1X Wash Buffer** as follows:
 - 9.7.1 Add 250 μ L of **1X Wash Buffer** to each assay well.
 - 9.7.2 Incubate for 30 seconds.
 - 9.7.3 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
 - 9.7.4 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
 - 9.7.5 Repeat steps 9.7.1 through 9.7.4 **four** more times.
- 9.8 Add 50 μ L of prepared **Detection Reagent A** and 50 μ L of **Detection Reagent B** to each well.
- 9.9 Gently mix the plate and incubate for 15 minutes **in the dark**.
- 9.10 Add 50 μ L of **Stop Solution** to each well and mix gently. Well color should change to yellow immediately. Add the **Stop Solution** in the same well order as done for the **Detection Reagents**.
- 9.11 Read the O.D. absorbance at 450 nm with a standard microplate reader within 5 minutes of stopping the reaction in step 9.10. If wavelength correction is available, set to 540 nm or 570 nm.

10. Calculation of Results

For analysis of the assay results, first derive the **Percent Relative Absorbance** for each sample or standard by calculating the mean absorbance between replicate wells. This value is then divided by the mean replicate Blank (0 ppb) absorbance, then multiplying by 100%:

$$\text{Percent Relative Absorbance (Sample or Standard)} = \frac{\text{Mean Absorbance (Sample or Standard)}}{\text{Mean Blank Absorbance}} \times 100\%$$

A standard curve is generated by plotting the **Percent Relative Absorbance** of each standard serial dilution point vs. the respective standard concentration. The Ciprofloxacin concentration contained in the samples can be interpolated by using linear regression of each sample **Percent Relative Absorbance** against the standard curve. This is best achieved using curve fitting software.


Note: If wavelength correction readings were available, subtract the readings at 540 nm or 570 nm from the readings at 450 nm. This may provide greater reading accuracy.

Note: If the samples measured were diluted, multiply the derived mean sample concentration by the dilution factor for a final sample concentration.

11. Typical Expected Data

11.1 Typical Data

For convenience in result calculation, absorbance as abscissa and standard conditions can be used as ordinate. The standard curve in the manual is only for reference; experimenters should draw the standard curve based on their own data.

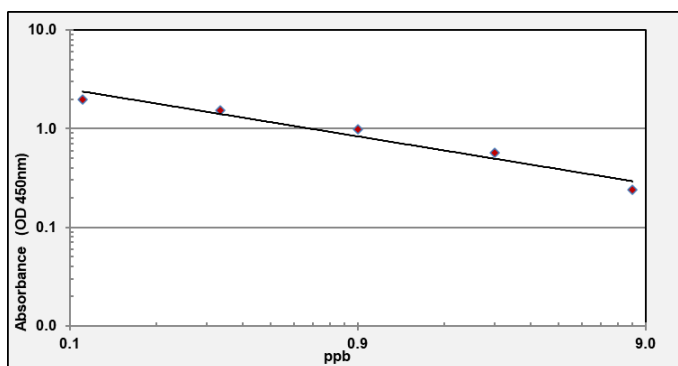

| ppb | Absorbance | | Mean Absorbance | B/B0(%) |
|-----|------------|-------|-----------------|---------|
| | Rep 1 | Rep 2 | | |
| 8.1 | 0.239 | 0.243 | 0.241 | 10.6 |
| 2.7 | 0.569 | 0.590 | 0.580 | 25.4 |
| 0.9 | 0.977 | 1.076 | 1.027 | 45.0 |
| 0.3 | 1.542 | 1.523 | 1.533 | 67.2 |
| 0.1 | 1.988 | 1.993 | 1.991 | 87.3 |
| 0.0 | 2.296 | 2.266 | 2.281 | 100.0 |

11.2 Limit of Detection

Derived by linear regression of OD450 of the Mean Blank + 2xSD

| Sample | Mean Recovery |
|-------------|---------------|
| Tissue | 0.3 ppb |
| Honey | 0.4 ppb |
| Milk | 3 ppb |
| Milk powder | 6 ppb |
| Egg | 3 ppb |

11.3 Reproducibility

Three samples of known concentrations were measured in 20 replicates within an assay plate to assess intra-assay reproducibility:

Mean Intra-Assay Precision - $\leq 8\%$ ($n = 3 \times 20$)

Three samples of known concentrations were measured across 40 replicate assays inter-assay reproducibility:

Mean Inter-Assay Precision - $\leq 10\%$ ($n = 3 \times 40$)

| Sample | Intra-Assay | | | Inter-Assay | | |
|------------|-------------|-------|-------|-------------|-------|-------|
| | 1 | 2 | 3 | 1 | 2 | 3 |
| Sample | 1 | 2 | 3 | 1 | 2 | 3 |
| n | 20 | 20 | 20 | 20 | 20 | 20 |
| Mean (ppb) | 0.335 | 1.079 | 1.701 | 0.297 | 1.281 | 1.626 |
| SD | 0.023 | 0.077 | 0.131 | 0.021 | 0.097 | 0.133 |
| CV (%) | 6.9 | 7.1 | 7.7 | 7.1 | 7.6 | 8.2 |

11.4 Cross-Reactivity Rate

| Substance | Cross Reactivity Rate |
|---------------|-----------------------|
| CPFX | 100% |
| Oxolinic Acid | 28% |
| Levofloxacin | 10% |
| Lomefloxacin | 4% |
| Marbofloxacin | 4% |
| Sarafloxacin | 2% |

11.5 Recovery

The recovery of Ciprofloxacin spiked at levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

| Sample | Number | Recovery Range |
|-------------|--------|----------------|
| Tissue | 10 | 85 \pm 15% |
| Honey | 10 | 85 \pm 10% |
| Milk | 10 | 85 \pm 15% |
| Milk powder | 10 | 85 \pm 15% |
| Egg | 10 | 85 \pm 10% |

12. Technical Resources

Technical Support:

For optimal service please be prepared to supply the lot number of the kit used.

USA

Aviva Systems Biology, Corp.
6370 Nancy Ridge Dr, Ste 104
San Diego, CA 92121

Phone: 858-552-6979
Toll Free: 888-880-0001
Fax: 858-552-6975

Technical support: techsupport@avivasysbio.com

China

Beijing AVIVA Systems Biology
6th Floor, B Building, Kaichi Tower
#A-2 Jinfu Road.
Daxing Industrial Development Zone
Beijing, 102600, CHINA

Phone: (86)10-60214720
Fax: (86)10-60214722
E-mail: support@avivasysbio.com.cn

中国地址: 北京大兴工业开发区金辅路甲 2 号凯驰大厦 B 座 6 层 (102600)

电话: 010-60214720/21

传真: 010-60214722

产品售前咨询及销售: sales@avivasysbio.com.cn

售后及技术支持: support@avivasysbio.com.cn