

SAA ELISA Kit (Mouse) (OKBB00384) Lot# KG0682

Instructions for Use

For the quantitative measurement of SAA in cell culture supernatants, serum, plasma (heparin, EDTA) and urine.

Variation between lots can occur. Refer to the manual provided with the kit.

This product is intended for research use only.



Table of Contents

1.	Background2
2.	Assay Summary3
3.	Storage and Stability3
4.	Kit Components3
5.	Precautions4
6.	Required Materials Not Supplied4
7.	Technical Application Tips4
8.	Reagent Preparation5
9.	Sample Preparation
10.	Assay Procedure8
11.	Calculation of Results9
12.	Typical Expected Data9
13.	Technical Resources



1. Background

Principle

Aviva Systems Biology SAA ELISA Kit (Mouse) (OKBB00384) is based on standard sandwich enzyme-linked immuno-sorbent assay technology. A mouse monoclonal antibody specific for SAA has been pre-coated onto 96-wellplate (12 x 8 Well Strips). Standards (E.coli; G20-Y122) and test samples are added to the wells and incubated. After washing, a biotinylated polyclonal goat detector antibody specific for SAA is added, incubated and followed by washing. Avidin-Biotin-Peroxidase Complex is then added, incubated and unbound conjugate is washed away. An enzymatic reaction is visualized through the addition of TMB substrate which is catalyzed by HRP to produce a blue color product that changes yellow after adding acidic stop solution. The density of yellow coloration read by absorbance at 450 nm and is quantitatively proportional to the amount of sample Mouse SAA captured in well.

Background


Serum amyloid A protein (SAA), also known as SAA1, is a protein that in humans is encoded by the SAA1 gene. This gene encodes a member of the serum amyloid A family of apolipoproteins. It is mapped to 11p15.1. SAA is a major acute phase protein that is highly expressed in response to inflammation and tissue injury. This protein also plays an important role in HDL metabolism and cholesterol homeostasis. High levels of this protein are associated with chronic inflammatory diseases including atherosclerosis, rheumatoid arthritis, Alzheimer's disease and Crohn's disease. This protein may also be a potential biomarker for certain tumors.

General Specifications

General Specifications						
Range 1.56 – 100 ng/mL						
Sensitivity < 0.15 ng/mL (Derived by linear regression of OD ₄₅₀ of the Mean Blank + 2xSD)						
	Mouse SAA					
	<u>UniProt ID</u> : P05366					
Specificity	<u>Gene ID</u> : 20208					
	Target Alias: Serum amyloid A-1 protein; Saa2; Saa-1					
Cross-Reactivity	No detectable cross-reactivity with other relevant proteins					



2. Assay Summary

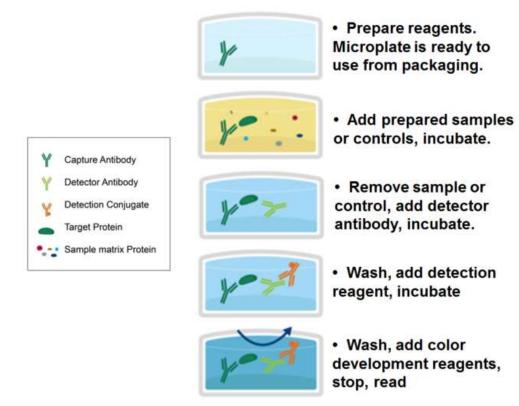

3. Storage and Stability

• Upon receipt store kit at -20°C until expiration date. Avoid multiple freeze/thaw cycles.

4. Kit Components

• The following reagents are the provided contents of the kit.

Description	Quantity	Storage Conditions
96-Well Plate Pre-coated with Anti-Mouse SAA Antibody	96 Wells (12 x 8 Well Strips)	
Lyophilized Recombinant Mouse SAA standard	2 x 100 ng	
100X Biotinylated Anti-Mouse SAA Antibody	1 x 100 μL	
100X Avidin-Biotin-Peroxidase Complex (ABC)	1 x 100 μL	
Sample Diluent Buffer	1 x 30 mL	Store at -20°C until
Antibody Diluent Buffer	1 x 12 mL	expiration date
ABC Diluent Buffer	1 x 12 mL	•
TMB Color Developing Agent	1 x 10 mL	
TMB Stop Solution	1 x 10 mL	
10X Wash Buffer	1 x 30 mL	



5. Precautions

- Read instructions fully prior to beginning use of the assay kit.
- Any deviations or modifications from the described method or use of other reagents could result in a reduction of performance.
- Reduce exposure to potentially harmful substances by wearing personal protective lab equipment including lab coats, gloves and glasses.
- For information on hazardous substances included in the kit please refer to the Material Safety Data Sheet (MSDS).

6. Required Materials Not Supplied

- Microplate reader capable of reading absorbance at 450 nm.
- · Automated plate washer (optional).
- Pipettes capable of precisely dispensing 0.5 µL through 1 mL volumes of aqueous solutions.
- Pipettes or volumetric glassware capable of precisely measuring 1 mL through 100 mL of aqueous solutions.
- New, clean tubes and/or micro-centrifuge tubes for the preparation of standards or samples.
- Absorbent paper or paper toweling.
- Distilled or deionized ultrapure water.

7. Technical Application Tips

- Do not mix or substitute components from other kits.
- To ensure the validity of experimental operation, it is recommended that pilot experiments using standards and a small selection of sample dilutions to ensure optimal dilution range for quantitation.
- Samples exhibiting OD measurements higher than the highest standard should be diluted further in the appropriate sample dilution buffers.
- Inspect all reagents prior to use. Components should contain no particulates or cloudiness and should be colorless.
- Prior to using the kit, briefly spin component tubes to collect all reagent at the bottom.
- Replicate wells are recommended for standards and samples.
- Cover microplate while incubating to prevent evaporation.
- Do not allow the microplate wells dry at any point during the assay procedure.
- Do not reuse tips or tube to prevent cross contamination.
- · Avoid causing bubbles or foaming when pipetting, mixing or reconstituting.
- Completely remove of all liquids when washing to prevent cross contamination.
- Prepare reagents immediately prior to use and do not store, with the exception of the top standard.
- Equilibrate all materials to ambient room temperature prior to use (standards exception).
- For optimal results in inter- intra-assay consistency, equilibrate all materials to room temperature prior to performing assay (standards exception) and perform all incubations at 37°C.
- Pipetting less than 1 µL is not recommended for optimal assay accuracy.
- Once the procedure has been started, all steps should be completed without interruption. Ensure that all reagents, materials and devices are ready at the appropriate time.
- Incubation times will affect results. All wells should be handled in the same sequential order and time intervals for optimal results.
- Samples containing precipitates or fibrin strands or which are hemolytic of lipemic might cause inaccurate results due to interfering factors.




8. Reagent Preparation

• Equilibrate all materials to room temperature prior to use and use immediately.

8.1 **SAA Assay Standards**

- 8.1.1 Prepare the SAA standards no greater than 2 hours prior to performing experiment. Standards should be held on ice until use in the experiment.
- 8.1.2 Reconstitute one of the provided 100 ng Lyophilized SAA Standard. Use one for each experiment. Prepare a stock 100 ng/mL SAA Standard by reconstituting one tube of Lyophilized SAA Standard as follows:
 - 8.1.2.1 Gently spin or tap the vial to collect all material at the bottom.
 - 8.1.2.2 Add 1 mL of **Sample Diluent Buffer** to the vial.
 - 8.1.2.3 Seal then mix gently and thoroughly.
 - 8.1.2.4 Leave the vial to sit at ambient temperature for 10 minutes.
- 8.1.3 Prepare a set of seven serially diluted standards as follows:
 - 8.1.3.1 Label tubes with numbers 1 8.
 - 8.1.3.2 Use the undiluted 100 ng/mL SAA Standard as the high standard point (Tube #1).
 - 8.1.3.3 Add 300 μ L of **Sample Diluent** to Tube #2 8.
 - 8.1.3.4 Prepare **Standard #2** by adding 300 μ L of 100 ng/mL **SAA** (Tube #1) to Tube #2. Mix gently and thoroughly.
 - 8.1.3.5 Prepare **Standard #3** by adding 300 μL of **Standard #2** from Tube #2 to Tube #3. Mix gently and thoroughly.
 - 8.1.3.6 Prepare further serial dilutions through Tube #7. Reference the table below as a guide for serial dilution scheme.
 - 8.1.3.7 Tube #8 is a blank standard (only **Sample Diluent Buffer**), which should be included with every experiment.

Standard Number (Tube)	Sample To Dilute	Volume Standard (μL)	Volume Sample Diluent Buffer (μL)	Total Volume (μL)	Final Concentration
1	100 ng/mL of SAA Standard	NA	1,000	1,000	100 ng/mL
2	100 ng/mL	300	300	600	50 ng/mL
3	50 ng/mL	300	300	600	25 ng/mL
4	25 ng/mL	300	300	600	12.5 ng/mL
5	12.5 ng/mL	300	300	600	6.25 ng/mL
6	6.25 ng/mL	300	300	600	3.12 ng/mL
7	3.12 ng/mL	300	300	600	1.56 ng/mL
8	NA	0	300	300	0.0 (Blank)

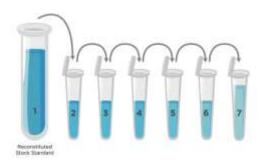



8.2 1X Biotinylated Anti-Mouse SAA Antibody

- 8.2.1 Prepare the **1X Biotinylated Anti-Mouse SAA Antibody** immediately prior to use by diluting the **100X Biotinylated Anti-Mouse SAA Antibody** 1:100 with Antibody Diluent Buffer.
- 8.2.2 For each 8 well strip to be used in the experiment prepare 1,000 μL by adding 10 μL of **100X** Biotinylated Anti-Mouse SAA Antibody to 990 μL Antibody Diluent Buffer.
- 8.2.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure.

8.3 1X Avidin-Biotin-Peroxidase Complex (ABC)

- 8.3.1 Prepare the **1X Avidin-Biotin-Peroxidase Complex (ABC)** immediately prior to use by diluting the **100X Avidin-Biotin-Peroxidase Complex (ABC)** 1:100 with **ABC Dilution Buffer**.
- 8.3.2 For each 8 well strip to be used in the experiment prepare 1,000 μL by adding 10 μL of **100X** Avidin-Biotin-Peroxidase Complex (ABC) to 990 μL ABC Dilution Buffer.
- 8.3.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure.

8.4 1X Wash Buffer

- 8.4.1 Add 270 mL of ultra-pure water to a clean > 500 mL bottle or other vessel.
- 8.4.2 Add the entire 30 mL contents of the **10X Wash Buffer** bottle to the water.
- 8.4.3 Seal and mix gently by inversion. Avoid foaming or bubbles.
- 8.4.4 Store the **1X Wash Buffer** at room temperature until ready to use in the procedure. Store the prepared **1X Wash Buffer** at 4°C for no longer than 1 week. Do not freeze.

8.5 Microplate Preparation

- Micro-plates are provided ready to use and do not require rinsing or blocking.
- Unused well strips should be returned to the original packaging, sealed and stored at 4°C.
- Equilibrate microplates to ambient temperatures prior to opening to reduce potential condensation.



9. Sample Preparation

9.1 Sample Preparation and Storage

- Store samples to be assayed at 2-8°C for 24 hours prior being assayed.
- For long term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.
- · Clear samples by centrifugation as follows:
 - Cell culture supernatant Remove particulates by centrifugation, analyze immediately or aliquot and store at -20°C
 - **Serum** Use a serum separator tube (SST) and allow serum to clot at room temperature for about four hours. Then, centrifuge for 15 min at approximately 1,000 x g. assay immediately or store samples at -20°C.
 - **Plasma** Collect plasma using heparin or EDTA as an anticoagulant. Centrifuge for 15 min at approximately 1,000 x g. Assay immediately or store samples at -20°C.
 - **Urine** Collect the first urine of the day, micturate directly into a sterile container. Remove impurities by centrifugation, assay immediately or aliquot and store samples at -20°C.

9.2 Sample Dilution

Target protein concentration must be estimated and appropriate sample dilution selected such that the final target protein concentration falls near the middle of the assay linear dynamic range.

- Prepare 150 µL sample for each replicate to be assayed.
- Dilute samples with Sample Diluent Buffer.
- Mix diluted samples gently and thoroughly.
- Pipetting less than 2 µL is not recommended for optimal assay accuracy.
- Refer to the following table for recommended sample dilution guidelines based on the dynamic range of this kit:

Estimated Sample T	arget Concentration	Dilution Level	Sample Volume For Two Replicates	Sample Diluent Buffer For Two Replicates	
High Concentration	100-1,000 ng/mL	1:100	3 µL	297 μL	
Medium Concentration	10-100 ng/mL	1:10	30 μL	270 μL	
Low Concentration	156-10,000 pg/mL	1:2	150 μL	150 μL	
Very Low Concentration	≤156 pg/mL	1:2 or No Dilution	-	-	



10. Assay Procedure

- The ABC Working Solution, TMB Color Developing Agent and TMB Stop Solution must be kept warm at 37°C for 30 minutes prior to use. Equilibrate all other reagents and materials to ambient room temperature prior to use in the procedure.
- To control for small potential variations in micro well-plate and day to day ambient temperature fluctuations, equilibrate all reagents prior to use and perform all incubation steps at to 37°C for optimal consistency and reproducibility.
- **10.1** Add 100 μL of serially titrated standards, diluted samples or blank into wells of the pre-coated well plate. At least two replicates of each standard, sample or blank is recommended.
- **10.2** Cover the plate with the plate sealer and incubate at 37°C for 90 minutes.
- **10.3** Remove the plate sealer and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
- **10.4** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 10.5 Add 100 µL of prepared 1X Biotinylated Anti-Mouse SAA Antibody to each well.
- **10.6** Cover with the plate sealer and incubate at 37°C for 60 minutes.
- **10.7** Wash plate 3 times with **1X Wash Buffer as** follows:
 - 10.7.1 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
 - 10.7.2 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
 - 10.7.3 Add 300 µL of **1X Wash Buffer** to each assay well.
 - 10.7.4 Incubate for 1 minute.
 - 10.7.5 Repeat steps 10.7.1 through 10.7.4 two more times.
- **10.8** Add 100 μL of prepared **1X Avidin-Biotin-Peroxidase Complex (ABC) into** each well, cover with plate sealer and incubate at 37°C for 30 minutes.
- **10.9** Wash plate 5 times with **1X Wash Buffer as** follows:
 - 10.9.1 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
 - 10.9.2 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
 - 10.9.3 Add 300 µL of 1X Wash Buffer to each assay well.
 - 10.9.4 Incubate for 1 minute.
 - 10.9.5 Repeat steps 10.9.1 through 10.9.4 four more times.
- **10.10** Add 90 μL of **TMB Color Developing Agent** to each well, cover with plate sealer and incubate at 37°C **in the dark** for 15-25 minutes.
 - (NOTE: optimal incubation time must be determined by the user. Optimal development can be visualized by blue shading in the top four standard wells, while the remaining standards are still clear.)
- 10.11 Add 100 µL of TMB Stop Solution to each well. Well color should change to yellow immediately.
- **10.12** Read the O.D. absorbance at 450 nm with a standard microplate reader within 30 minutes of stopping the reaction in step 10.11.



11. Calculation of Results

For analysis of the assay results, calculate the **Relative OD**₄₅₀ for each test or standard well as follows:

(Relative
$$OD_{450}$$
) = (Well OD_{450}) – (Mean Blank Well OD_{450})

The standard curve is generated by plotting the mean replicate **Relative OD**₄₅₀ of each standard serial dilution point vs. the respective standard concentration. The **SAA** concentration contained in the samples can be interpolated by using linear regression of each mean sample **Relative OD**₄₅₀ against the standard curve. This is best achieved using curve fitting software.


Note: if the samples measured were diluted, multiply the derived mean sample concentration by the dilution factor for a final sample concentration.

12. Typical Expected Data

12.1 Typical absorbance values. Expected absorbance for standards when TMB incubation is performed for 20 minutes at 37°C and measured at OD₄₅₀.

Standard Number	8	7	6	5	4	3	2	1
Concentration (ng/mL)	0	1.56	3.12	6.25	12.5	25	50	100
OD ₄₅₀	0.079	0.108	0.215	0.437	0.535	0.990	1.423	1.872

12.2 <u>Typical standard curve</u>. This standard curve is for demonstration purposes only. An assay specific standard curve should be performed with each assay.

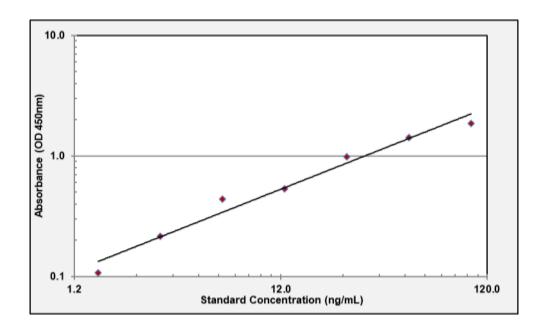



12.3 Reproducibility

		Intra-Assay		Inter-Assay			
Sample ID	1	2	3	1	2	3	
n =	16	16	16	24	24	24	
Mean Measured Concentration (pg/mL)	3290	11528	49554	3575	12179	49101	
Standard Deviation (pg/mL)	151.34	691.68	2775.02	178.75	937.78	2749.65	
Consistency (%CV)	4.6	6	5.6	5	7.7	5.6	



13. Technical Resources

13.1 References

- 13.1.1 "Entrez Gene: SAA1 serum amyloid A1".
- 13.1.2 Glaser T, Housman D, Lewis WH, Gerhard D, Jones C (Jan 1990). "A fine-structure deletion map of human chromosome 11p: analysis of J1 series hybrids". Somat Cell Mol Genet 15 (6): 477–501.
- 13.1.3 Polyak K, Xia Y, Zweier JL, Kinzler KW, Vogelstein B (Sep 1997). "A model for p53-induced apoptosis". Nature 389 (6648): 300–5.

13.2 Technical Support

<u>USA</u>

Aviva Systems Biology, Corp. 10211 Pacific Mesa Blvd, Ste 401 San Diego, CA 92121

Phone: 858-552-6979 Toll Free: 888-880-0001 Fax: 858-552-6975

Technical support: techsupport@avivasysbio.com

China

Beijing AVIVA Systems Biology 6th Floor, B Building, Kaichi Tower #A-2 Jinfu Road. Daxing Industrial Development Zone Beijing, 102600, CHINA

Phone: (86)10-60214720 Fax: (86)10-60214722

E-mail: support@avivasysbio.com.cn

中国地址:北京大兴工业开发区金辅路甲2号凯驰大厦B座6层(102600)

电话: 010-60214720/21 传真: 010-60214722

产品售前咨询及销售: sales@avivasysbio.com.cn

售后及技术支持: support@avivasysbio.com.cn