



**Papillomavirus L1-Capsids Antibody (IgG)  
(Human) ELISA Kit  
(OKCA00473)  
Instructions for use**

For the quantitative measurement of Papillomavirus L1-Capsids Antibody (IgG) in serum and plasma.

This product is intended for research use only.

Lot to lot kit variations can occur. Use the kit manual which has been provided along with the kit packaging.

## Table of Contents

1. Background .....	2
2. Assay Summary .....	2
3. Precautions .....	4
4. Storage and Stability .....	3
5. Kit Components .....	3
6. Required Materials Not Supplied .....	4
7. Technical Application Tips .....	4
8. Reagent Preparation .....	5
9. Sample Preparation .....	7
10. Assay Procedure .....	8
11. Calculation of Results .....	9
12. Typical Expected Data .....	9
13. Technical Resources .....	11

## 1. Background


### Principle

Aviva Systems Biology Papillomavirus L1-Capsids Antibody (IgG) (Human) ELISA Kit (OKCA00473) is based on standard reverse capture sandwich enzyme-linked immuno-sorbent assay technology. Papillomavirus L1-Capsids antigen has been pre-coated and blocked in a 96-wellplate (12 x 8 Well Strips). Standards or test samples are added to the wells, incubated and washed. An HRP conjugated detector antibody specific for Human IgG is added, incubated and followed by washing. An enzymatic reaction is produced through the addition of substrate which is catalyzed by HRP generating a blue color product that changes to yellow after adding acidic stop solution. The density of yellow coloration is read by absorbance at 450 nm and is quantitatively proportional to the amount of sample Papillomavirus L1-Capsids Antibody (IgG) captured in well.

### General Specifications

General Specifications	
Specificity	Human Papillomavirus L1-Capsids Antibody (IgG)
Cross-Reactivity	No detectable cross-reactivity with other relevant proteins

## 2. Assay Summary

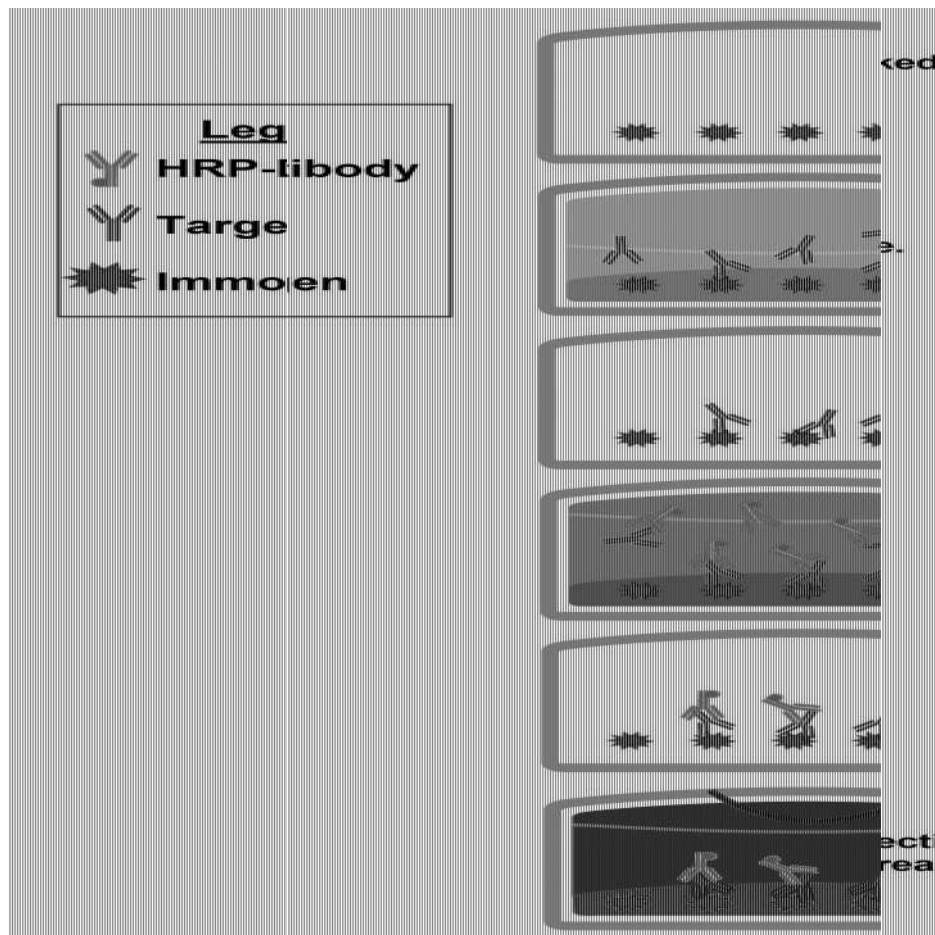

### 3. Storage and Stability

- Upon receipt store kit at 4°C for 6 months.

### 4. Kit Components

- The following reagents are the provided contents of the kit.

Description	Quantity	Storage Conditions
Papillomavirus L1-Capsids Antibody (IgG) Microplate	96 Wells (12 x 8 Well strips)	4°C for 6 Months
Negative Control	1 x 800 µl	
Positive Control	1 x 800 µl	
100X Anti-Human IgG HRP-conjugate	1 x 120 µl	
HRP-conjugate Diluent	1 x 20 ml	
Sample Diluent	2 x 20 ml	
25X Wash Buffer	1 x 20 ml	
TMB Substrate	1 x 10 ml	
Stop Solution	1 x 10 ml	

### 5. Precautions

- Read instructions fully prior to beginning use of the assay kit.
- Any deviations or modifications from the described method or use of other reagents could result in a reduction of performance.
- Reduce exposure to potentially harmful substances by wearing personal protective lab equipment including lab coats, gloves and glasses.
- For information on hazardous substances included in the kit please refer to the Material Safety Data Sheet (MSDS).
- Kit cannot be used beyond the expiration date on the label.

### 6. Required Materials Not Supplied

- Microplate reader capable of reading absorbance at 450 nm.
- Automated plate washer (optional).
- Pipettes capable of precisely dispensing 0.5 µL through 1 mL volumes of aqueous solutions.
- Pipettes or volumetric glassware capable of precisely measuring 1 mL through 100 mL of aqueous solutions.
- New, clean tubes and/or micro-centrifuge tubes for the preparation of standards or samples.
- Absorbent paper or paper toweling.
- Distilled or deionized ultrapure water.
- 37°C Incubator (optional)

## 7. Technical Application Tips

- Do not mix or substitute components from other kits.
- To ensure the validity of experimental operation, it is recommended that pilot experiments using standards and a small selection of sample dilutions to ensure optimal dilution range for quantitation.
- Samples exhibiting OD measurements higher than the highest standard should be diluted further in the appropriate sample dilution buffers.
- Prior to using the kit, briefly spin component tubes to collect all reagents at the bottom.
- Replicate wells are recommended for standards and samples.
- Cover microplate while incubating to prevent evaporation.
- Do not allow the microplate wells dry at any point during the assay procedure.
- Do not reuse tips or tube to prevent cross contamination.
- Avoid causing bubbles or foaming when pipetting, mixing or reconstituting.
- Completely remove of all liquids when washing to prevent cross contamination.
- Prepare reagents immediately prior to use and do not store, with the exception of the top standard.
- Equilibrate all materials to ambient room temperature prior to use (standards exception).
- For optimal results for inter- and intra-assay consistency, equilibrate all materials to 37°C prior to performing assay (standards exception) and perform all incubations at 37°C.
- Pipetting less than 1  $\mu$ L is not recommended for optimal assay accuracy.
- Once the procedure has been started, all steps should be completed without interruption. Ensure that all reagents, materials and devices are ready at the appropriate time.
- Incubation times will affect results. All wells should be handled in the same sequential order and time intervals for optimal results.
- Samples containing precipitates, fibrin strands or bilirubin, or are hemolytic or lipemic might cause inaccurate results due to interfering factors.
- TMB Substrate is easily contaminated and should be colorless or light blue until added to plate. Handle carefully and protect from light.

## 8. Reagent Preparation

- Equilibrate all materials to room temperature prior to use and use prepare immediately prior to use.
- The following reagents are provided at ready to use concentrations and require no preparation:
  - Negative Control
  - Positive Control
  - Anti-Human IgG HRP-Conjugate
  - HRP-Conjugate diluent
  - Sample Diluent
  - TMB Substrate
  - Stop Solution

### 8.1 1X Anti-Human IgG HRP-Conjugate

- 8.1.1 Prepare the **1X Anti-Human IgG HRP-Conjugate** immediately prior to use by diluting the **100X Anti-Human IgG HRP-Conjugate** 1:100 with **HRP-Conjugate Diluent**.
- 8.1.2 For each well strip to be used in the experiment (8-wells) prepare 1,000  $\mu\text{L}$  of **1X Anti-Human IgG HRP-Conjugate** by adding 10  $\mu\text{L}$  of **100X Anti-Human IgG HRP-Conjugate** to 990  $\mu\text{L}$  **HRP-Conjugate Diluent**.
- 8.1.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure. Do not store at 1X concentration for future use.

### 8.2 1X Wash Buffer

- 8.2.1 If crystals have formed in the 25X **Wash Buffer** concentrate, equilibrate to room temperature and mix gently until crystals have completely dissolved.
- 8.2.2 Add the entire 20 mL contents of the 25X **Wash Buffer** bottle to 480 mL of ultra-pure water to a clean > 1,000 mL bottle or other vessel.
- 8.2.3 Seal and mix gently by inversion. Avoid foaming or bubbles.
- 8.2.4 Store the **1X Wash Buffer** at room temperature until ready to use in the procedure. Store the prepared **1X Wash Buffer** at 4°C for no longer than 1 week. Do not freeze.

### 8.3 Microplate Preparation

- Micro-plates are provided ready to use and do not require rinsing or blocking.
- Unused well strips should be returned to the original packaging, sealed and stored at 4°C.
- Equilibrate microplates to ambient temperatures prior to opening to reduce potential condensation.

## 9. Sample Preparation

### 9.1 Sample Preparation and Storage

- Store samples to be assayed at 2-8°C for 24 hours prior being assayed.
- For long term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.
- Samples not indicated in the manual must be tested to determine if the kit is valid.
- Prepare samples as follows:
- Serum - Use a serum separator tube (SST) and allow samples to clot for two hours at room temperature or overnight at 4°C before centrifugation for 15 minutes at 1,000 x g. Remove serum and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- Plasma - Collect plasma using EDTA, or heparin as an anticoagulant. Centrifuge for 15 minutes at 1,000 x g at 2-8°C within 30 minutes of collection. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.

### 9.2 Sample Dilution (1:1000)

Target protein concentration must be estimated and appropriate sample dilution selected such that the final target protein concentration falls near the middle of the assay linear dynamic range. Samples exhibiting saturation should be further diluted.

- Dilute samples using **Sample Diluent**.
- Mix diluted samples gently and thoroughly.
- Pipetting less than 2 µL is not recommended for optimal assay accuracy.

Serum and plasma samples require a 1000-fold dilution into Sample Diluent. The suggested 1000-fold dilution can be achieved by adding 5µl sample to 95µl of Sample Diluent first, then complete the 1000-fold dilution by adding 5µl of this solution to 245µl of Sample Diluent.

## 10. Assay Procedure

- Equilibrate all reagents and materials to ambient room temperature prior to use in the procedure.
  - Optimal results for inter- and intra-assay reproducibility will be observed when incubation at 37°C as indicated in the procedure below.
- 10.1** Determine the required number of wells and return any remaining unused wells and desiccant to the pouch.
  - 10.2** Assign a well for an absolute blank.
  - 10.3** Add 100 µL of the **1:1000 diluted samples, Positive Control** and **Negative Control** into wells of the **Papillomavirus L1-Capsids Antibody (IgG) Microplate**. At least two replicates of each control, sample or blank is recommended.
  - 10.4** Cover the plate with the well plate lid and incubate at 37°C for 30 minutes.
  - 10.5** Remove the plate lid and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
  - 10.6** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
  - 10.7** Wash plate 3 times with **1X Wash Buffer** as follows:
    - 10.7.1 Add 300 µL of **1X Wash Buffer** to each assay well.
    - 10.7.2 Incubate for 2 minute.
    - 10.7.3 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
    - 10.7.4 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
    - 10.7.5 Repeat steps 10.7.1 through 10.7.4 **two** more times.
  - 10.8** Add 100 µL of the **1X Anti-Human IgG HRP-Conjugate** to all the wells.
  - 10.9** Cover with the well-plate lid and incubate at 37°C for 30 minutes.
  - 10.10** Remove the plate lid and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
  - 10.11** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
  - 10.12** Repeat the wash in step 10.7 with **five** washes.
  - 10.13** Add 90 µL of **TMB Substrate** to each well and incubate at 37°C **in the dark** for 20 minutes. Wells should change to gradations of blue. If the color is too deep, reduce the incubation time.  
(NOTE: optimal incubation time must be determined by the user. Optimal development can be visualized by blue shading in the top four standard wells, while the remaining standards are still clear.)
  - 10.13** Add 50 µL of **Stop Solution** to each well. Well color should change to yellow immediately. Add the **Stop Solution** in the same well order as done in step 10.13. Gently tap the plate to mix.
  - 10.14** Read the O.D. absorbance at 450 nm with a standard microplate reader within 10 minutes of stopping the reaction in step 10.14. If wavelength correction is available, set to 540 nm or 570 nm.



## 11 Calculation of Results

A positive or negative Papillomavirus L1-Capsids Antibody (IgG) determination is derived by comparing the test samples to the **Positive** and **Negative Controls**.

While  $OD_{\text{sample}} < 2.1 \times OD_{\text{negative}}$ : Negative

While  $OD_{\text{sample}} \geq 2.1 \times OD_{\text{negative}}$ : Positive

## 12 Typical Expected Data

### 12.13 Reproducibility

Intra-assay Precision: 3 samples with known low, middle and high levels Papillomavirus L1-Capsids Antibody (IgG) were tested with 20 replicates on one plate, respectively.

Inter-assay Precision: 3 samples with known low, middle and high level Papillomavirus L1-Capsids Antibody (IgG) were tested on 3 different plates, 8 replicates in each plate.

Mean Intra-Assay:  $CV \leq 15\%$

Mean Inter-Assay:  $CV \leq 15\%$

## 13 Technical Resources

### Technical Support:

For optimal service please be prepared to supply the lot number of the kit used.

#### USA

Aviva Systems Biology, Corp.  
5754 Pacific Center Blvd, Suite 201  
San Diego, CA 92121

Phone: 858-552-6979  
Toll Free: 888-880-0001  
Fax: 858-552-6975

Technical support: [techsupport@avivasysbio.com](mailto:techsupport@avivasysbio.com)

#### China

Beijing AVIVA Systems Biology  
6th Floor, B Building, Kaichi Tower  
#A-2 Jinfu Road.  
Daxing Industrial Development Zone  
Beijing, 102600, CHINA

Phone: (86)10-60214720  
Fax: (86)10-60214722  
E-mail: [support@avivasysbio.com.cn](mailto:support@avivasysbio.com.cn)

中国地址: 北京大兴工业开发区金辅路甲 2 号凯驰大厦 B 座 6 层 (102600)  
电话: 010-60214720/21  
传真: 010-60214722

产品售前咨询及销售: [sales@avivasysbio.com.cn](mailto:sales@avivasysbio.com.cn)  
售后及技术支持: [support@avivasysbio.com.cn](mailto:support@avivasysbio.com.cn)