

# LGALS3 ELISA Kit (Dog) (OKEH00478)

# **Instructions for Use**

For the quantitative measurement of LGALS3 in serum, plasma, tissue homogenates, cell culture supernatants and other biological fluids.

This product is intended for research use only.



# Table of Contents

| 1.  | Background                      | . 2 |
|-----|---------------------------------|-----|
| 2.  | Assay Summary                   | . 2 |
| 3.  | Precautions                     | . 4 |
| 4.  | Storage and Stability           | . 3 |
| 5.  | Kit Components                  | . 3 |
| 6.  | Required Materials Not Supplied | . 4 |
| 7.  | Technical Application Tips      | . 4 |
| 8.  | Reagent Preparation             | . 5 |
| 9.  | Sample Preparation              | . 7 |
| 10. | Assay Procedure                 | . 8 |
| 11. | Calculation of Results          | . 9 |
| 12. | Typical Expected Data           | . 9 |
| 13. | Technical Resources             | 10  |



## 1. Background

#### **Principle**

Aviva Systems Biology LGALS3 ELISA Kit (Dog) (OKEH00478) is based on standard sandwich enzyme-linked immuno-sorbent assay technology. An antibody specific for LGALS3 has been pre-coated onto a 96-well plate (12 x 8 Well Strips). Standards or test samples are added to the wells, incubated and removed. A biotinylated detector antibody specific for LGALS3 is added, incubated and followed by washing. Avidin-Peroxidase Conjugate is then added, incubated and unbound conjugate is washed away. An enzymatic reaction is produced through the addition of TMB substrate which is catalyzed by HRP generating a blue color product that changes to yellow after adding acidic stop solution. The density of yellow coloration read by absorbance at 450 nm and is quantitatively proportional to the amount of sample LGALS3 captured in the well.

#### Target Background


Galactose-specific lectin which binds IgE. May mediate with the alpha-3, beta-1 integrin the stimulation by CSPG4 of endothelial cells migration. Together with DMBT1, required for terminal differentiation of columnar epithelial cells during early embryogenesis. In the nucleus: acts as a pre-mRNA splicing factor. Involved in acute inflammatory responses including neutrophil activation and adhesion, chemoattraction of monocytes macrophages, opsonization of apoptotic neutrophils, and activation of mast cells.

#### **General Specifications**

| General Specifications |  |  |  |  |  |  |  |
|------------------------|--|--|--|--|--|--|--|
| Range                  | 78 - 5000 pg/mL  |  |  |  |  |  |  |
| LOD                    | NA   |  |  |  |  |  |  |
| Specificity            | Dog LGALS3<br>UniProt ID: P38486<br>GeneID: 404021<br>Target Alias: Lectin L-29, Laminin-binding protein, IgE-binding protein, Galactose-specific lectin 3,<br>Mac-2 antigen, 35 kDa lectin, Carbohydrate-binding protein 35(CBP 35) |  |  |  |  |  |  |
| Cross-Reactivity       | No detectable cross-reactivity with other relevant proteins  |  |  |  |  |  |  |



## 2. Assay Summary

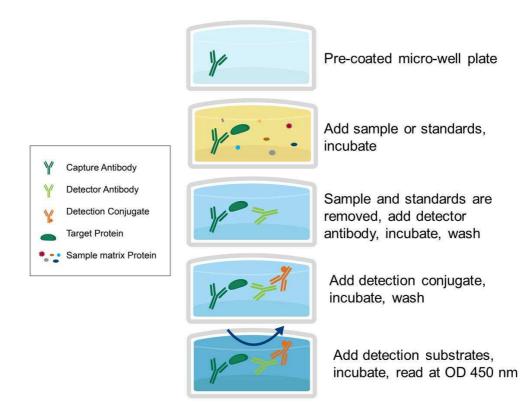

## 3. Storage and Stability

• Upon receipt store kit at 4°C for 6 months or -20°C for long term storage (exceptions noted below). Avoid any freeze/thaw cycles.

## 4. Kit Components

• The following reagents are the provided contents of the kit.

| Description                                | Quantity                      | Storage Conditions  |  |
|--|-------------------------------|---------------------|--|
| LGALS3 Microplate                          | 96 Wells (12 x 8 Well strips) |                     |  |
| LGALS3 Lyophilized Standard                | 2                             |                     |  |
| 100X Biotinylated LGALS3 Detector Antibody | 1 x 120 μL                    |                     |  |
| 100X Avidin-HRP Conjugate                  | 1 x 120 μL                    | 4°C for 6 Months    |  |
| Sample Diluent                             | 1 x 20 mL                     | -20°C for 12 Months |  |
| Detector Antibody Diluent                  | 1 x 10 mL                     |                     |  |
| Conjugate Diluent                          | 1 x 10 mL                     |                     |  |
| 25X Wash Buffer                            | 1 x 30 mL                     |                     |  |
| Stop Solution                              | 1 x 10 mL                     | Store at 4°C        |  |
| TMB Substrate                              | 1 x 10 mL                     |                     |  |



## 5. Precautions

- Read instructions fully prior to beginning use of the assay kit.
- Any deviations or modifications from the described method or use of other reagents could result in a reduction of performance.
- Reduce exposure to potentially harmful substances by wearing personal protective lab equipment including lab coats, gloves and glasses.
- For information on hazardous substances included in the kit please refer to the Material Safety Data Sheet (MSDS).
- Kit cannot be used beyond the expiration date on the label.

## 6. Required Materials Not Supplied

- Microplate reader capable of reading absorbance at 450 nm.
- Automated plate washer (optional).
- Pipettes capable of precisely dispensing 0.5 µL through 1 mL volumes of aqueous solutions.
- Pipettes or volumetric glassware capable of precisely measuring 1 mL through 100 mL of aqueous solutions.
- New, clean tubes and/or micro-centrifuge tubes for the preparation of standards or samples.
- Absorbent paper or paper toweling.
- Distilled or deionized ultrapure water.
- 37°C Incubator (optional)

## 7. Technical Application Tips

- Do not mix or substitute components from other kits.
- To ensure the validity of experimental operation, it is recommended that pilot experiments using standards and a small selection of sample dilutions to ensure optimal dilution range for quantitation.
- Samples exhibiting OD measurements higher than the highest standard should be diluted further in the appropriate sample dilution buffers.
- Prior to using the kit, briefly spin component tubes to collect all reagents at the bottom.
- Replicate wells are recommended for standards and samples.
- Cover microplate while incubating to prevent evaporation.
- Do not allow the microplate wells dry at any point during the assay procedure.
- Do not reuse tips or tube to prevent cross contamination.
- Avoid causing bubbles or foaming when pipetting, mixing or reconstituting.
- Completely remove of all liquids when washing to prevent cross contamination.
- Prepare reagents immediately prior to use and do not store, with the exception of the top standard.
- Equilibrate all materials to ambient room temperature prior to use (standards exception).
- For optimal results for inter- and intra-assay consistency, equilibrate all materials to 37°C prior to performing assay (standards exception) and perform all incubations at 37°C.
- Pipetting less than 1 µL is not recommended for optimal assay accuracy.
- Once the procedure has been started, all steps should be completed without interruption. Ensure that all reagents, materials and devices are ready at the appropriate time.
- Incubation times will affect results. All wells should be handled in the same sequential order and time intervals for optimal results.
- Samples containing bilirubin, precipitates or fibrin strands or are hemolytic or lipemic might cause inaccurate results due to interfering factors.
- TMB Substrate is easily contaminated and should be colorless or light blue until added to plate. Handle carefully and protect from light.




## 8. Reagent Preparation

• Equilibrate all materials to room temperature prior to use and use prepare immediately prior to use.

#### 8.1 Dog LGALS3 Assay Standards

- 8.1.1 Prepare the LGALS3 standards no greater than 2 hours prior to performing experiment. Standards should be held on ice until use in the experiment.
- 8.1.2 Reconstitute one vial of the provided 5 ng Lyophilized LGALS3 Standard for each experiment. Prepare the stock 5,000 pg/mL Standard by reconstituting one tube of Lyophilized LGALS3 Standard as follows:
  - 8.1.2.1 Gently spin or tap the vial at 6,000 10,000 rpm for 30 seconds to collect all material at the bottom.
  - 8.1.2.2 Add 1 mL of Sample Diluent to the vial.
  - 8.1.2.3 Seal the vial then mix gently and thoroughly.
  - 8.1.2.4 Leave the vial at ambient temperature for 15 minutes.
- 8.1.3 Prepare a set of seven serially diluted standards as follows:
  - 8.1.3.1 Label tubes with numbers 2 8.
  - 8.1.3.2 Use the undiluted 5,000 pg/mL LGALS3 Standard as the high standard point (Tube #1).
  - 8.1.3.3 Add 300  $\mu L$  of Sample Diluent to Tube #2 8.
  - 8.1.3.4 Prepare **Standard #2** by adding 300 μL of 5,000 pg/mL **LGALS3** (Tube #1) to Tube #2. Mix gently and thoroughly.
  - 8.1.3.5 Prepare **Standard #3** by adding 300 μL of **Standard #2** from Tube #2 to Tube #3. Mix gently and thoroughly.
  - 8.1.3.6 Prepare further serial dilutions through Tube #7. Reference the table below as a guide for serial dilution scheme.
  - 8.1.3.7 Tube #8 is a blank standard (only **Sample Diluent**), which should be included with every experiment.

| Standard<br>Number<br>(Tube) | Standard To Dilute                              | Volume Standard<br>to Dilute<br>(µL) | Volume Sample<br>Diluent Buffer<br>(μL) | Total Volume<br>(μL) | Final Concentration |
|------------------------------|---|--------------------------------------|---|----------------------|---------------------|
| 1                            | 5,000 pg/mL<br>Reconstituted<br>LGALS3 Standard | 1,000                                | NA                                      | NA                   | 5,000 pg/mL         |
| 2                            | 5,000 pg/mL                                     | 300                                  | 300                                     | 600                  | 2,500 pg/mL         |
| 3                            | 2,500 pg/mL                                     | 300                                  | 300                                     | 600                  | 1,250 pg/mL         |
| 4                            | 1,250 pg/mL                                     | 300                                  | 300                                     | 600                  | 625 pg/mL           |
| 5                            | 625 pg/mL                                       | 300                                  | 300                                     | 600                  | 312 pg/mL           |
| 6                            | 312 pg/mL                                       | 300                                  | 300                                     | 600                  | 156 pg/mL           |
| 7                            | 156 pg/mL                                       | 300                                  | 300                                     | 600                  | 78 pg/mL            |
| 8                            | NA  | 0                                    | 300                                     | 300                  | 0.0 (Blank)         |

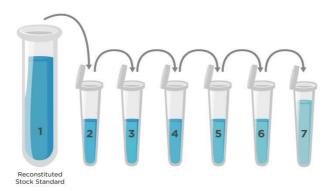



#### 8.2 1X Biotinylated LGALS3 Detector Antibody

- 8.2.1 Prepare the **1X Biotinylated LGALS3 Detector Antibody** immediately prior to use by diluting the **100X Biotinylated LGALS3 Detector Antibody** 1:100 with **Detector Antibody Diluent**.
- 8.2.2 For each well strip to be used in the experiment (8-wells) prepare 800 μL by adding 8 μL of 100X Biotinylated LGALS3 Detector Antibody to 792 μL Detector Antibody Diluent. Prepare approximately 200 μL extra for pipetting overage each time the assay is performed.
- 8.2.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure. Do not store at 1X concentration for future use.

#### 8.3 1X HRP-Avidin Conjugate

- 8.3.1 Prepare the **1X Avidin-HRP Conjugate** immediately prior to use by diluting the **100X Avidin-HRP Conjugate** 1:100 with **Conjugate Diluent**.
- 8.3.2 For each well strip to be used in the experiment (8-wells) prepare 800 μL by adding 8 μL of 100X Avidin-HRP Conjugate to 792 μL Conjugate Diluent. Prepare approximately 200 μL extra for pipetting overage each time the assay is performed.
- 8.3.3 Mix thoroughly and gently. Hold no longer than 2 hours prior to using in procedure. Do not store at 1X concentration for future use.

#### 8.4 1X Wash Buffer

- 8.4.1 If crystals have formed in the **25X Wash Buffer** concentrate, equilibrate to room temperature and mix gently until crystals have completely dissolved.
- 8.4.2 Add the entire 30 mL contents of the **25X Wash Buffer** bottle to 720 mL of ultra-pure water to a clean > 1,000 mL bottle or other vessel.
- 8.4.3 Seal and mix gently by inversion. Avoid foaming or bubbles.
- 8.4.4 Store the **1X Wash Buffer** at room temperature until ready to use in the procedure. Store the prepared **1X Wash Buffer** at 4°C for no longer than 1 week. Do not freeze.

#### 8.5 Microplate Preparation

- Micro-plates are provided ready to use and do not require rinsing or blocking.
- Unused well strips should be returned to the original packaging, sealed and stored at °4C.
- Equilibrate microplates to ambient temperatures prior to opening to reduce potential condensation.



## 9. Sample Preparation

#### 9.1 Sample Preparation and Storage

- Store samples to be assayed at 2-8°C for 24 hours prior being assayed.
- For long term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.
- Samples not indicated in the manual must be tested to determine if the kit is valid.
- Prepare samples as follows:
  - Serum Use a serum separator tube (SST) and allow samples to clot for 30 minutes at room temperature or overnight at 4°C before centrifugation for 15 minutes at 1,000 x g. Remove serum and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
  - **Plasma** Collect plasma using EDTA, or heparin as an anticoagulant. Centrifuge for 15 minutes at 1,000 x g at 2-8°C within 30 minutes of collection. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
  - Tissue Homogenates Rinse 100 mg of tissue with 1X PBS then homogenize in 20 ml of 1X PBS and store overnight at -20°C. Perform two freeze-thaw cycles to break the cell membranes, then centrifuge for 5 minutes at 5,000 x g, 2-8°C. Remove the supernatant and assay immediately. Alternatively, aliquot and store samples at -20°C or -80°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.
  - Cell culture supernatants and other biological fluids Remove particulates by centrifugation and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze/thaw cycles.

#### 9.2 Sample Dilution

Target protein concentration must be estimated and appropriate sample dilution selected such that the final target protein concentration falls near the middle of the assay linear dynamic range. Samples exhibiting saturation should be further diluted.

- Dilute samples using Sample Diluent.
- Mix diluted samples gently and thoroughly.
- Pipetting less than 2 µL is not recommended for optimal assay accuracy.



## **10. Assay Procedure**

- Equilibrate all reagents and materials to ambient room temperature prior to use in the procedure.
- Optimal results for intra- and inter-assay reproducibility will be obtained when performing incubation steps at 37°C as indicated below.
- **10.1** Determine the required number of wells and return any remaining unused wells and desiccant to the pouch.
- **10.2** Add 100 µL of serially titrated standards, diluted samples or blank into wells of the LGALS3 Microplate. At least two replicates of each standard, sample or blank is recommended.
- **10.3** Cover the plate with the well plate lid and incubate at 37°C for 2 hours.
- **10.4** Remove the plate lid and discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
- **10.5** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 10.6 Add 100 µL of prepared 1X Biotinylated LGALS3 Detector Antibody to each well.
- 10.7 Cover with the well-plate lid and incubate at 37°C for 60 minutes.
- **10.8** Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
- **10.9** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- 10.10 Wash plate 3 times with 1X Wash Buffer as follows:
  - 10.10.1 Add 300 µL of **1X Wash Buffer** to each assay well.
  - 10.10.2 Incubate for 1 minute.
  - 10.10.3 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle.
  - 10.10.4 Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
  - 10.10.5 Repeat steps 10.10.1 through 10.10.4 **two** more times.
- **10.11** Add 100 μL of prepared **1X Avidin-HRP Conjugate** into each well and incubate at 37°C for 60 minutes.
- 10.12 Discard the liquid in the wells by rigorously flicking into an acceptable waste receptacle or aspiration.
- **10.13** Gently blot any remaining liquid from the wells by tapping inverted on the benchtop onto paper toweling. Do not allow the wells to completely dry at any time.
- **10.14** Wash plate **5 times** with **1X Wash Buffer** as in Step 10.10.
- **10.15** Add 90 μL of **TMB Substrate** to each well and incubate at 37°C in the dark for 15-30 minutes. Wells should change to gradations of blue. If the color is too deep, reduce the incubation time.

(NOTE: optimal incubation time must be determined by the user. Optimal development can be visualized by blue shading in the top four standard wells, while the remaining standards are still clear.)

- **10.16** Add 50 μL of **Stop Solution** to each well. Well color should change to yellow immediately. Add the **Stop Solution** in the same well order as done for the **TMB Substrate**.
- **10.17** Read the O.D. absorbance at 450 nm with a standard microplate reader within 5 minutes of stopping the reaction in step 10.16. If wavelength correction is available, set to 540 nm or 570 nm.



### **11. Calculation of Results**

For analysis of the assay results, calculate the Relative OD<sub>450</sub> for each test or standard well as follows:

(Relative  $OD_{450}$ ) = (Well  $OD_{450}$ ) – (Mean Blank Well  $OD_{450}$ )

The standard curve is generated by plotting the mean replicate **Relative OD**<sub>450</sub> of each standard serial dilution point vs. the respective standard concentration. The LGALS3 concentration contained in the samples can be interpolated by using linear regression of each mean sample **Relative OD**<sub>450</sub> against the standard curve. This is best achieved using curve fitting software.

**Note:** if wavelength correction readings were available, subtract the readings at 540 nm or 570 nm from the readings at 450 nm. This may provide greater reading accuracy.

**Note:** if the samples measured were diluted, multiply the derived mean sample concentration by the dilution factor for a final sample concentration.

## 12. Typical Expected Data


#### 12.1 Reproducibility

Three samples concentrations were measured in replicate within an assay plate and across replicate assays to assess Intra- and Mean Inter-Assay Precision.

Mean Intra-Assay Precision - ≤ 4.9%, n = 20

Mean Inter-Assay Precision - ≤ 7.5%, n = 20

**12.2 Typical standard curve**. This standard curve is for demonstration purposes only. An assay specific standard curve should be performed with each assay.

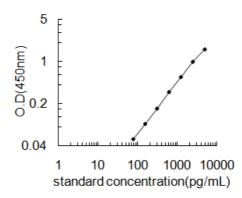

#### 12.3 Recovery

The recovery of LGALS3 spiked at levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Mean recovery when spiking into Serum and Plasma = 101%



## **13. Technical Resources**

#### **Technical Support:**

For optimal service please be prepared to supply the lot number of the kit used.

#### <u>USA</u>

Aviva Systems Biology, Corp. 5754 Pacific Center Blvd, Suite 201 San Diego, CA 92121

Phone: 858-552-6979 Toll Free: 888-880-0001 Fax: 858-552-6975

Technical support: techsupport@avivasysbio.com

<u>China</u>

Beijing AVIVA Systems Biology 6th Floor, B Building, Kaichi Tower #A-2 Jinfu Road. Daxing Industrial Development Zone Beijing, 102600, CHINA

Phone: (86)10-60214720 Fax: (86)10-60214722 E-mail: support@avivasysbio.com.cn

中国地址:北京大兴工业开发区金辅路甲2号凯驰大厦B座6层(102600) 电话:010-60214720/21 传真:010-60214722

产品售前咨询及销售: sales@avivasysbio.com.cn 售后及技术支持: support@avivasysbio.com.cn